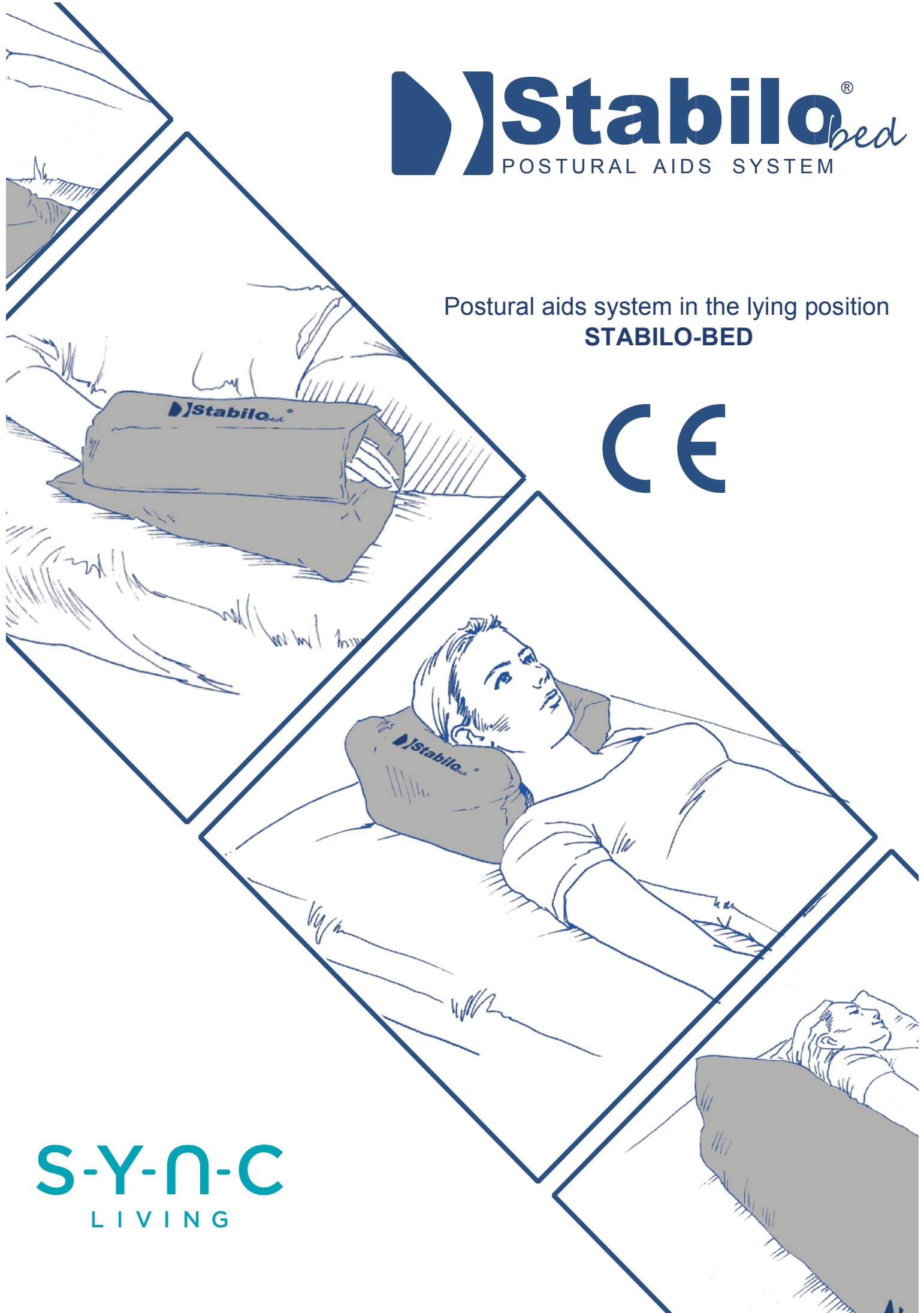


Stabilo[®]bed

POSTURAL AIDS SYSTEM

Postural aids system in the lying position
STABILO-BED

CE








S-Y-n-C
LIVING







STABILO

Version: 16.05.2018

STABILObed® is the modular 1st medical class system that is used for precise and effective positioning the patient in the lying position. **STABILObed®** system is a range of blocks and cushions, used for stabilization of different body parts. They are covered by innovative fabric **PureMed™**. This material is very safe for human use and easy to sterilize.

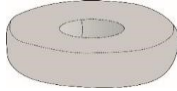





















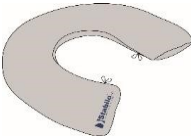
FILLING



















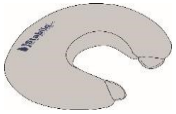
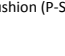
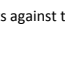
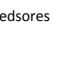







	<p>Granules – micro-beads with small granulometry. The micro-beads may go smoothly inside the cushion and because of that, it can be anatomically shaped very easily. What is more, the PureMed™ cover, with high vapour permeability, allows the air to transfer through the cushion.</p>
	<p>VISmemo™ - is our innovative memory foam. This is the most active in areas with high bed sores risk.</p>
	<p>Comfort foam – soft foam with a lot of closed, small air bubbles to maintain the patient's shape.</p>
	<p>Off-loading foam – rigid technical foam with a closed structure allowing the patient to return to their initial position even after long-term lying.</p>
	<p>Silicone – this filling is very soft and lightweight and can maintain the perfect anatomic shape. It's Pressure-resistant even after long-term lying (supine) pressure.</p>





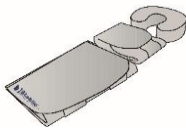
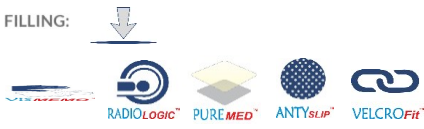

	<p>PureMed™ is very safe and easy to sterilize. It consists of antibacterial components based on silver zeolite (AgZ) and because of it, there is a possibility of multiple use by patients in hospitals or rehabilitation centers.</p>
	<p>Each component of STABILObed® is covered by soft, skin-friendly material.</p>
	<p>The components can be used to immobilise patients during diagnostic imaging (MRI, CT scan, X-ray).</p>
	<p>Some of the STABILObed® elements are covered by special non slip fabric, preventing against unwanted moving elements on the mattress.</p>
	<p>STABILObed® components are equipped with Velcro straps to prevent their movement when placed on a self-adhesive sheet.</p>
	<p>Certain components contain a special system of adjustment located in the cover which helps to shape the STABILObed® parts.</p>

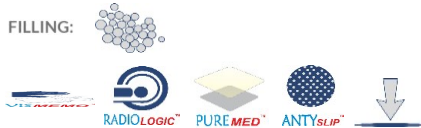

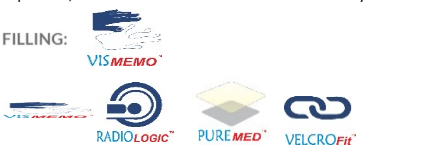

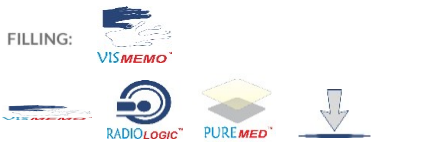
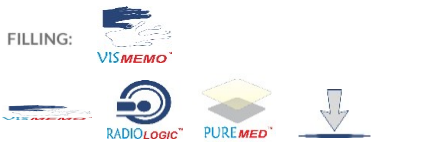
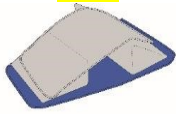
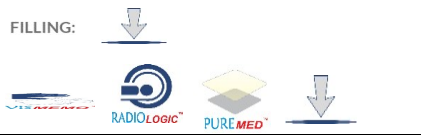

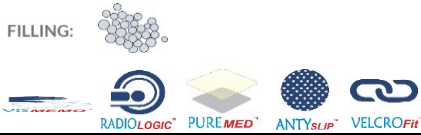
STABILObed®




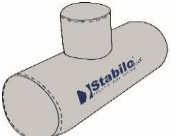


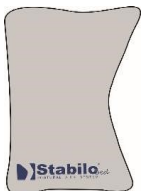
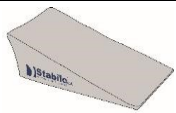


PRODUCT	DESCRIPTION	CODE
<p>P-SS-01 LOWER EXTREMITY ABDUCTION PILLOW</p> 	<ul style="list-style-type: none"> Our abduction pillow protects the patient from adduction contracture of the hips and, furthermore, it decreases the risk of bedsores that can appear on medial side of lower limbs. <p>FILLING: </p> <p>FEATURES:    </p>	<p>P-SS-01/1</p>
<p>P-SS-02 CYLINDRICAL CUSHION</p> <p>P-SS-02</p> 	<ul style="list-style-type: none"> Our bolster knee roll pillow can be used in various ways, but used under knees allows to achieve Semi- Fowler position. Cylindrical cushion should be used by patients who are placed in bed for a long time. It brings the lower back into a zero-gravity position, improves circulation and prevents against the swollen lower extremity. Also, when you use our abduction block (P-SS-03) you will maintain the legs in anti-bedsores position (in abduction). <p>FILLING: </p> <p>FEATURES:  </p>	<p>P-SS-02/1</p> <p>P-SS-02/2</p>
<p>P-SS-06 CALF WEDGE</p> <p>P-SS-06</p> 	<ul style="list-style-type: none"> It prevents from bedsores that can appear on heel bone. Reduces risk of thrombophlebitis (vein inflammation). The length is excellent to support the knee against the hyperextension. <p>FILLING: </p> <p>FEATURES:    </p>	<p>P-SS-06/1</p>
<p>P-SS-08 LATERAL POSITIONING CUSHION</p> <p>P-SS-08</p> 	<ul style="list-style-type: none"> Keeps the patient in the lateral (safe) position and supports knees, pelvic and upper body. Allows patient to keep semi-ventral and semi-lateral position. Lateral positioning cushion can be used with thermoactive support cushion (P-SS-22), which prevents against the bedsores additionally (look at pages: 35-36). Highly recommended for pregnant women to maintain safe position on the left side for long night. Available in one size compatible with all patient sizes. There is special string and ribbing inside that allow to achieve a required position of the cushion. <p>FILLING: </p> <p>FEATURES:   </p>	<p>P-SS-08/1</p>

<p>P-SS-09 ANTI BEDSORE ROUND CUSHION P-SS-09</p>	<ul style="list-style-type: none"> • Our anti-bedsore round cushion is in decubitus prophylaxis and in bedsores, postoperative and after injury treatment. • When it's used under head, it prevents against the bedsores on the occipital bone and support the cervical lordosis. • The inner hole in the cushion effectively relieves the sensitive points of the patient's body. It doesn't allow them to come into contact with a hospital mattress, what improves the blood supply and nourishments of soft tissues. 	<p>P-SS-09/1</p>
	<p>FILLING: </p> <p>FEATURES:    </p>	<p>P-SS-09/2</p> <p>P-SS-09/3</p>
<p>P-SS-10 KNEE WEDGE P-SS-10</p>	<ul style="list-style-type: none"> • Using of the knee wedge in laying position reduces the pressure and pain in the lower back. • Used together with ring cushion (P-SS-13) achieves Semi-Fowler position, where pelvis and lower back is at 30°. • It improves blood circulation and reduces the risk of lower limb swelling. • Prevents against the knee hyperextension. <p>FILLING: </p> <p>FEATURES:    </p>	<p>P-SS-10/1</p>
<p>P-SS-12 LATERAL POSITIONING WEDGE P-SS-12</p>	<ul style="list-style-type: none"> • Highly recommended for patients who are placed in bed for a long time and required frequent changes of position. • It provides support for the trunk from one side. • Relieves the sensitive areas of the body and significantly reduces the risk of bedsores. • Available in one size compatible with all patient sizes. 	<p>P-SS-12/1</p>
	<p>FILLING: </p> <p>FEATURES:    </p>	<p>P-SS-12/2</p>
<p>P-SS-13 RING CUSHION P-SS-13</p>	<ul style="list-style-type: none"> • Long roller is ideal for patients requiring body stabilization in lying position. • It helps in proper head positioning and contributes to anti-bedsore prophylaxis by relieving sensitive areas of the body. • It can also be used by pregnant women to increase the comfort of lying down. • It is also possible to purchase thermoactive support cushion. • There is special string and ribbing inside that allow to achieve a required position of the cushion. <p>FILLING: </p> <p>FEATURES:    </p>	<p>P-SS-13/1</p>
		
<p>P-SS-14 SEMI-SOFT BOLSTER PILLOWS P-SS-14</p>	<ul style="list-style-type: none"> • The semi-soft bolster pillows are the perfect solution to place the patient in semi-lateral or semi-ventral positions. <ul style="list-style-type: none"> • They improve anti-bedsore prophylaxis. • Using of semi-soft bolster pillows helps to prevent against the limbs deformation and muscles contractures. • To increase anti-bedsore prophylaxis it is highly recommend to use bolster pillow together with thermoactive 	<p>P-SS-14/1</p>

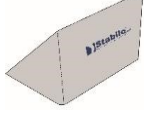

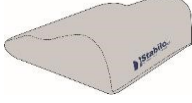

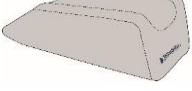

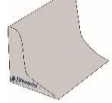
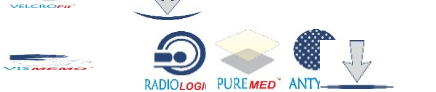
	<p>support cushion (P-SS-22).</p> <p>FILLING:  FEATURES:  </p>	<p>P-SS-14/2</p>
<p>P-SS-15 V PILLOW P-SS-15</p>	<ul style="list-style-type: none"> • Our 'V' shaped cushion comfortably supports the head and prevents against the bedsores. • It is great aid in sitting and semi-reclined position. • Our 'V' pillow can be used with thermoactive support cushion. 	<p>P-SS-15/1</p>
	<p>FILLING:  FEATURES:  </p>	
<p>P-SS-17 SEMI-RIGID BOLSTER PILLOWS FOR TRUNK SUPPORT P-SS-17</p>	<ul style="list-style-type: none"> • The semi-rigid bolster pillows are the perfect solution to place the patient in lateral, semi-lateral or semi-ventral positions. <ul style="list-style-type: none"> • They improve anti-bedsores prophylaxis. • We highly recommend to use our semi-rigid bolster pillow with shaped wedges (K-SS-08), that prevent against the slipping and changing position of bolster pillow. 	<p>P-SS-17/1</p>
	<p>FILLING:  FEATURES:  </p>	<p>P-SS-17/2</p>
<p>P-SS-18 HEADREST P-SS-18</p>	<ul style="list-style-type: none"> • Headrest makes it possible to place the lying patient in middle position. • It prevents from cervical spinal deformations and prevents against the bedsores. 	<p>P-SS-18/1</p>
	<p>FILLING:  FEATURES:    </p>	
<p>P-SS-19 NECK RING CUSHION P-SS-19</p>	<ul style="list-style-type: none"> • Our neck ring cushion allows to place the head in the middle position both in sitting and semi-reclined positions. • It supports the cervical spine and prevents against the neck deformations. • It prevents against the bedsores on the occipital bone. • Neck ring cushion can be used with thermoactive support cushion (P-SS-22), with prevents against the bedsores additionally 	<p>P-SS-19/1</p>
	<p>FILLING:  FEATURES:  </p>	
<p>P-SS-20 UPPER EXTREMITY SUPPORT WITH PROTECTOR P-SS-20</p>	<ul style="list-style-type: none"> • It's equipped with original memory foam VISmemo™, that prevents patient from bedsores that can appear in elbow and forearm area. • It stabilizes upper limb and helps to avoid the muscle contractions. • Available in one size compatible with all patient sizes. <p>FILLING:  </p>	<p>P-SS-20/1</p>
	<p>FEATURES:    </p>	

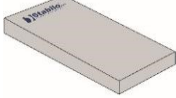

<p>P-SS-21 SILICONE PILLOW P-SS-21</p>	<ul style="list-style-type: none"> The silicone pillow helps to put the patient in a safe position and improves anti-bedsore prophylaxis. It prevents against the pathological deformities and neck muscles contractures. The cushion is covered with a soft, skin-friendly material in white colour. To place the patient in lateral or semi-lateral lying position, there is necessary to use the self-adhesive wedges (K-SS-02) or shaped wedges (K-SS-08) together with silicone pillow 	<p>P-SS-21/1</p>
	<p>FILLING:</p> 	
<p>P-SS-22 THERMOACTIVE SUPPORT CUSHION P-SS-22</p>	<ul style="list-style-type: none"> Thermoactive support cushion is made of memory foam VISmemo™. This foam moulds to the body in response to heat and pressure, evenly distributing body weight. It then returns to its original shape once you remove the pressure. This high solution increases a blood and oxygen supply in soft tissues. It's highly recommended for decubitus prophylaxis and it should be used in areas of the body at high risk for pressure Sores. Thermoactive support cushion is available with others parts of STABLObed® system. 	<p>P-SS-22/1</p>
	<p>FILLING:</p> 	
<p>P-SS-23 PEDIATRIC IMMOBILIZER P-SS-23</p>	<ul style="list-style-type: none"> Immobilizer is designed for children with spasticity. It secures children with limbs and trunk deformities. It prevents deepening deformities and helps in the process of physiotherapy. It can also be used to transport a child with recent trauma. The product consists of three elements: the head stabilizer and 2 independent elements with Velcro straps. Immobilizer has no steel elements so it can be used in X- ray, MRI and CT scanning procedures. 	<p>P-SS-23/1</p>
	<p>FILLING:</p> 	<p>P-SS-23/2</p>
<p>P-SS-24 ANTI-BEDSORE DONUT CUSHION P-SS-24</p>	<ul style="list-style-type: none"> This anti-bedsore cushion is used in the anti-bedsore prophylaxis and treatment of the existing bedsores as well as in post-operative and post-traumatic conditions. When used on the seat of a stroller or a car, it protects the patient against the bedsores in the area of the ischial tubers and provides an adequate distribution of body weight. The inner hole in the cushion effectively relieves the sensitive points of the patient's body, which are properly nourished and supplied with blood when do not have a direct contact with the seat. 	<p>P-SS-24/1</p>
	<ul style="list-style-type: none"> In order to maintain a proper hygiene of the product, it is recommended to use an additional terry cover. 	<p>P-SS-24/2</p>

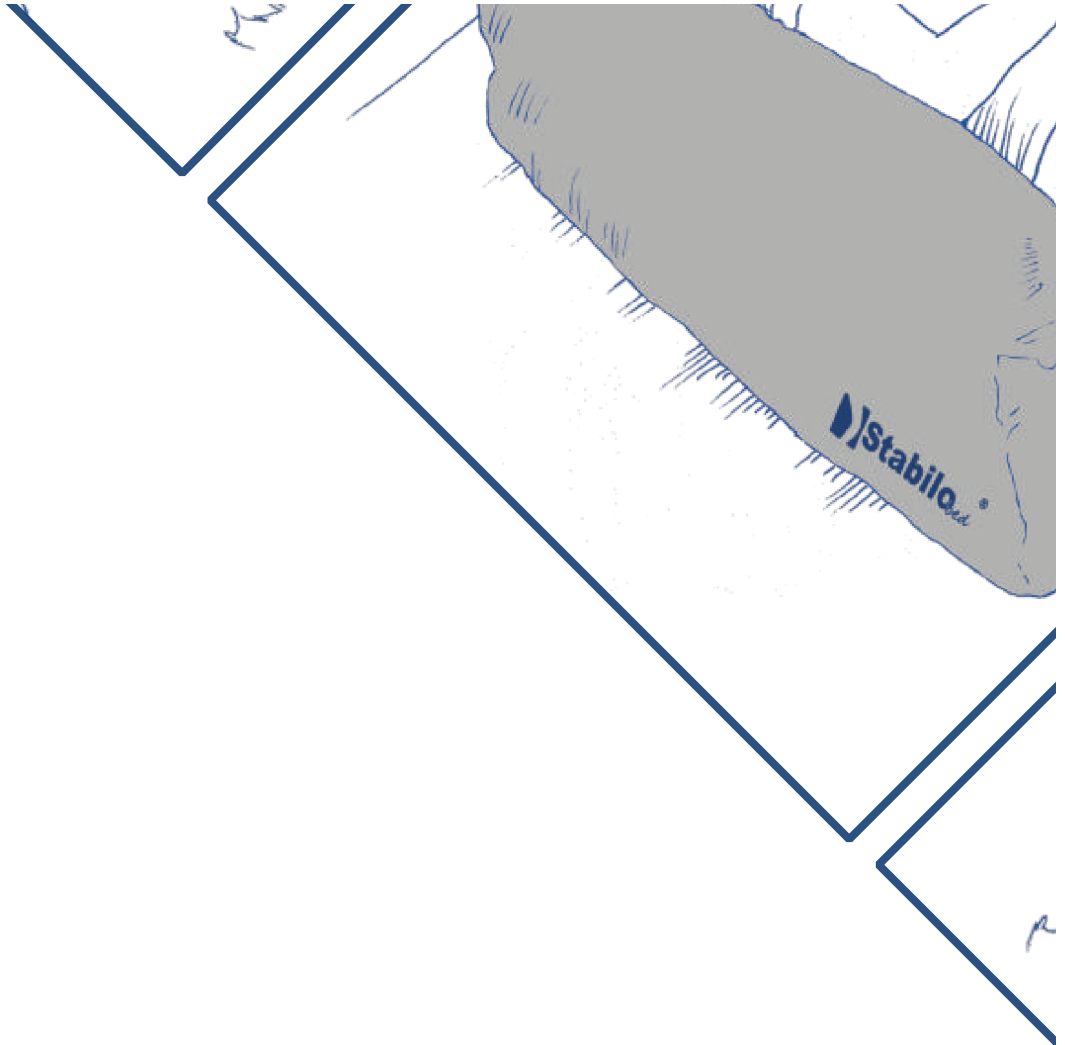
	<p>FILLING:</p> 	<p>P-SS-24/3</p>
<p>P-SS-25 HEAD CUSHION P-SS-25</p> 	<p>Thanks to its design, the cushion provides an additional side support for the head, and the filling with the thermoplastic shape memory foam VISMEMO™ is a guarantee of the highest comfort while lying down. The cushion reduces the pressure in the occiput area which significantly decreases the risk of bedsores. Also, it is equipped with a special Velcro strip to prevent the cushion from slipping during using. In order to keep the proper hygiene of the product, it is recommended to use an additional terry cover.</p> <p>FILLING:</p> 	<p>P-SS-25/1</p>
<p>P-SS-26 ORTHOPEDIC NECK PILLOW P-SS-26</p> 	<p>This pillow has been designed and made as an anatomically profiled 3D pillow. It is made of memory shape foam VISMEMO™ which adjusts perfectly to the patient's body shape. Thanks to its design, it provides a proper support for the neck and head. The thermoelastic shape memory foam keeps the shape of the body for a short time after we stop putting pressure on it and slowly returns to its original form. The cushion is equipped with a cover made of special PureMed™ fabric, which is very easy to disinfect.</p> <p>FILLING:</p> 	<p>P-SS-26/1</p>
	<p>FILLING:</p> 	<p>P-SS-26/2</p>
<p>P-SS-27 REFLUX BABY WEDGE P-SS-27</p> 	<ul style="list-style-type: none"> The anti-reflux wedge has been created for infants with problems related to the gastroesophageal reflux. Using the wedge significantly relieves the effects of the condition by putting the little patients in the proper lateral position and securing them at the same time. The correct position of the child prevents the backflow of the gastroesophageal content, which significantly reduces the risk of choking. <p>FILLING:</p> 	<p>P-SS-27/1</p>
<p>P-SS-28 PEDIATRIC LOWER LIMBS ABDUCTION WEDGE P-SS-28</p> 	<ul style="list-style-type: none"> Paediatric lower limbs abduction wedge protects the young patients against the adduction contractures of the hip joints and reduces the risk of bedsores on the medial side of the lower limbs. It helps to maintain the correct supine position by controlling the position of the legs. In order to prevent leg crossing, the wedge has been equipped with a special upper cover fastened with VelcroFit™ <p>FILLING:</p> 	<p>P-SS-28/1</p>

<p>P-SS-29 UNIVERSAL POSITIONING CUSHION P-SS-29</p> 	<ul style="list-style-type: none"> The cushion is a perfect addition in post-operative therapy. Its size and appropriate hardness helps in proper support of different body parts, e.g. a foot in numerous therapeutic positions. The cover is made of patient-friendly and easy to disinfect PureMed™ material and AntySlip™ material, also used in the cover, prevents the cushion from moving. <p>FILLING: </p> 	<p>P-SS-29/1</p>
<p>P-SS-30 LOWER LIMBS POSITIONING WEDGE P-SS-30</p>	<ul style="list-style-type: none"> Thanks to its design, the wedge enables the correct positioning of the lower limbs in the bends, at the same protecting the patient from the adduction contractures in the hip joints. Using the cushion prevents the knee hyperextension in the knee joints and reduces the risk of bedsores on the inner parts of the legs. It is highly recommended for the bedridden patients. It can be successfully used in both the supine and semi-reclining position, e.g. in post-operative therapy. 	<p>P-SS-30/1</p>
	<p>FILLING: </p> 	<p>P-SS-30/2</p>
<p>P-SS-31 UNIVERSAL CUSHION P-SS-31</p> 	<ul style="list-style-type: none"> A cushion with an anatomic neck cut stabilizes the head in the proper position as well as fills the cervical lordosis and prevents the bedsores around the occiput area. The versatility of the cushion makes it possible to keep the patient in a neutral position and reduces the spinal pain. When used between the legs, it protects the patient against the bedsores around knees and ankles. Using the cushion also reduces the risk of a pathological adduction, e.g. after the endo prosthetic hip replacement. The cushion has anti-bedsores and stabilizing functions. It can be used in numerous ways, such as: <ul style="list-style-type: none"> - an under head cushion, - a cushion used between or under the lower limbs, a support for upper limb elevation (anti-swelling and anti bedsores function), - an aid in maintaining the semi-lateral and semi recumbent position on the stomach. In order to maintain a proper hygiene of the product, it is recommended to use an additional terry cover. 	<p>P-SS-31/1</p>
<p>K-SS-01 THERMOACTIVE ELEVATION ARM WEDGE K-SS-01</p>	<ul style="list-style-type: none"> Elevation arm wedge stabilizes the arm in rising position and decrease the swelling of hand and forearm. It's equipped with original memory foam VISmemo™, that prevents patient from bedsores that can appear in elbow and forearm area. It's highly recommended for women after mastectomy or other surgeries where axillary lymph nodes were excised. Available in one size compatible with all patient sizes. 	<p>K-SS-01/1</p>
	<p>FILLING: </p> 	



<p>K-SS-02 WEDGE WITH SELF-GRIPPING FUNCTION K-SS-02</p>	<ul style="list-style-type: none"> The wedge can be attach to the sheet, because of self-gripping function, and prevents against the slipping of cushion. Highly recommended for heavy users, patients with spasticity, as well as patients that remain in the supine position for a long time. It provides additional stabilization and comfort. 	<p>K-SS-02/1</p>
	<p>FILLING:</p> 	<p>K-SS-02/2</p>
<p>K-SS-06 THERMOACTIVE LOWER LIMB RAISING WEDGE K-SS-06</p>	<ul style="list-style-type: none"> The wedge is made of original memory foam VISmemo™. It maintains the lower limbs slightly raised, what improves blood circulation and reduces the risk of swollen ankles. It can be used as a support for compression therapy and anti-decubitus treatment. We recommend use it at night. Because of our raising wedge, the blood flows from legs to the heart, the lower extremity circumference is reduced what allows to wake up in the bed in the morning and putting the compression stockings on. Our thermoactive wedge is anatomic-shaped what prevents against the lower limb deformations. It maintains pelvic tilt and provides zero gravity back pain relief. 	<p>K-SS-06/1</p>
		<p>K-SS-06/2</p>
<p>K-SS-07 THERMOACTIVE BRAUN SPLINT K-SS-07</p>	<ul style="list-style-type: none"> Our Braun splint is made of original memory foam VISmemo™. It maintains the lower limb slightly raised, what improves blood circulation and reduces the risk of swollen ankle. It can be used as a support for compression therapy and anti-decubitus treatment. It prevents from bedsores that can appear on heel bone. Reduces risk of thrombophlebitis (vein inflammation). The length is excellent to support the knee against the hyperextension. Our thermoactive splint is anatomic-shaped what prevents against the lower limb deformations. 	<p>K-SS-07/1</p>
	<p>FILLING:</p> 	<p>K-SS-07/2</p>
<p>K-SS-08 SHAPED WEDGE K-SS-08</p>	<ul style="list-style-type: none"> It protects the semi-rigid bolster pillow(P-SS-17) from slipping. Especially recommended for heavy users and patients with spasticity. It provides additional stabilization and comfort. 	<p>K-SS-08/1</p>
		<p>K-SS-08/2</p>

<p>B-SS SELF-ADHESIVE SHEET B-SS</p>	<ul style="list-style-type: none"> • The sheet is made of soft terry in a grey or blue colour. • Due to self-adhesive function, it is possible to attach STABLObed® cushion to it in order to prevent its moving. • The sheet ends with gum surface to prevent slipping out of the mattress. 	<p>B-SS/1</p>
	<ul style="list-style-type: none"> • The sheet size matches the standard mattresses, there is also possibility of custom size sheets. 	<p>B-SS/2</p>
		<p>B-SS/3</p>



S-Y-N-C

LIVING

Unit 3, 1 York Road
Belfast
BT15 3GU

T: +4428 90-74 5333
E: sales@syncliving.co.uk
W: syncliving.co.uk